



ROTHBAND

PATIENT POSITIONING

INSTRUCTIONS FOR USE

LATERAL FOOT IMAGING STEP

VERSION 2.0

**UK
CA**

Table of Contents

INTRODUCTION	3
INTENDED USE	3
POSITIONING RECOMMENDATIONS	3
TRANSPORTATION	3
TECHNICAL INFORMATION	4
COMPLIANCE	4
ASSEMBLY INSTRUCTIONS	5
Parts List.....	5
Attaching the Handrail	5
CAUTIONS AND WARNINGS	6
CLEANING INSTRUCTIONS	7
CONTACT DETAILS	7

INTRODUCTION

These Instructions cover the following models:

PRODUCT

LFIS – Lateral Imaging Step

INTENDED USE

The Rothband Imaging Step range has been designed for use during x-ray examinations of the lower limb.

As all Radiographers are aware, there are many variations in technique both departmentally and due to clinical situation and patient build. These recommendations are intended as a guide to enable a local protocol to be established and may be refined and adapted as and when required.

POSITIONING RECOMMENDATIONS



Standing Lateral Foot



AP Ankle

A closed cell foam Lateral foot block is provided with each step. This may be attached to the step if desired, with the Velcro tape provided. Replacement blocks are available from Rothband.

TRANSPORTATION

The step is provided with wheels to facilitate its movement between examinations. Simply lift the end opposite the wheels until the step rides on the wheels for easy movement.

TECHNICAL INFORMATION

PRODUCT SPECIFICATIONS	
Product Dimensions	900 mm x 510 mm x 1540 mm
Product Weight	15 Kg
Weight Limit	150 Kg
Height to First Step	210 mm
Height from First Step to Second Step	200 mm
Dimensions of FB40L	360 mm x 30 mm x 100 mm (15mm Slot)

COMPLIANCE

The Lateral Foot Imaging Step complies with the Medical Device Directive (MDD) 93/42/EEC. and has been designed and manufactured to comply with the following standards:

Standard	Description
BS EN 14183:2003	Step Stools
EN ISO 14971:2012	Medical devices - Application of risk management to medical devices (ISO 14971:2007, Corrected version 2007-10-01)
EN ISO 15223-1:2016	Medical devices - Symbols to be used with medical device labels, labelling and information to be supplied - Part 1: General requirements.

ASSEMBLY INSTRUCTIONS


Parts List

- 12 x Allen Bolts with standard washers
- 1 x Allen Key
- 1 x Foam Block (FP40L)
- 1 x Base Unit
- 2 x Handle
- 1 x Cross Piece
- 2 x 8" Strips of Velcro
- 1 x Foam Care Guide

Attaching the Handrail

Mount each side handle using the bolts and shakeproof washers provided. There are four bolts per handle. An Allen key is provided for the bolts.

Fix the cross piece using 4 Allen bolts and shakeproof washers.

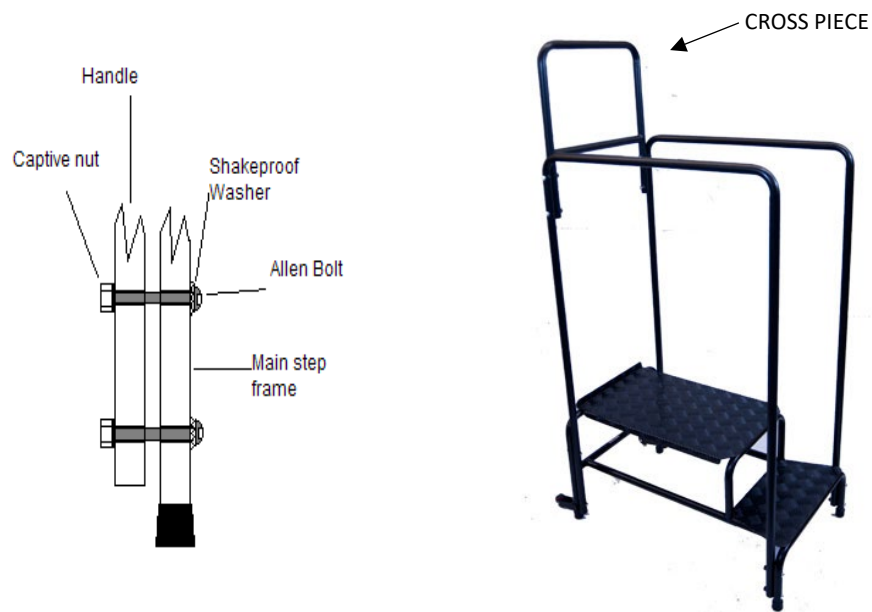



CAUTIONS AND WARNINGS

- Intended for professional healthcare facility environments where operators with medical training are continually available when patients are present.
- Individuals embarking the step should be closely supervised and assisted.
- Individuals should not be left un-supervised whilst on the step.
- There are no special skills required to use this positioning aid.
- The step is only suitable for indoor use.
- The step should only be used on a flat and level surface.
- The medical professional is responsible for applying his/her best medical judgment when positioning a patient.
- The step is suitable for continuous use.
- Not suitable for sterilisation.
- Do not place device on or near a heat source.
- Do not use abrasive cleaners, please see relevant section for the correct cleaning instructions.
- The step should be regularly inspected for damage and loose fasteners.
- The above warnings, cautions and any safety considerations should be observed on a routine and regular basis, not only upon installation.
- The device has a weight limit of 150Kg

CLEANING INSTRUCTIONS

Use a solution of mild soap or detergent or disinfectant* in hand-hot water.

Apply using a soft cloth & light pressure to prevent swirls.

Rinse off with clean water and dry by blotting with a soft cloth/towel.

WSR Medical Solutions Limited recommend Envirolex hard surface cleaner.

*DO NOT Use Solvents or Alcohol based cleaning products.

DO NOT Use any scouring materials/pads to clean the product as this will damage the surface finish.



Cross contamination is a major hazard to patients. The Imaging Steps should therefore be cleaned in between patients.

CONTACT DETAILS



WSR Medical Solutions
(Trading as Rothband)
Unit E Bridgewater Close,
Network 65 Business Park,
Burnley,
Lancashire,
BB11 5TE,
United Kingdom

Tel: 01706 830086

e-mail: sales@rothband.com